
DIGITAL DISPLAY ALCOHOL BREATH TESTER

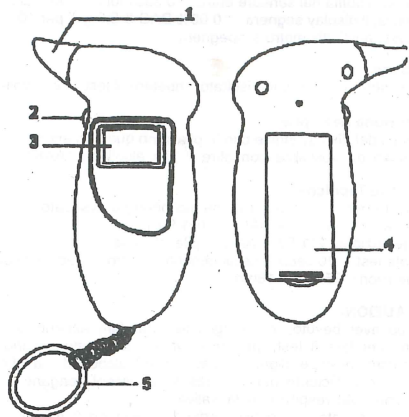
- (1) Breath inhaler
(2) Power button
(3) LCD display
(4) Battery compartment

Features

1. Digital alcohol test with audible alert
2. Range: 0.00-0.19% BAC & 0.0-1.9g/l
3. Quick response and resume
4. Auto power off
5. Use 2x "AAA" alkaline batteries

Installing battery

1. Slide open the battery compartment lid from the back of the unit.

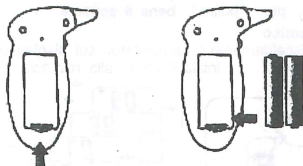
2. Insert two "AAA" (UM-4) alkaline batteries with correct polarity as indicated.
3. Replace the battery compartment lid.

How to use the alcohol tester

Warm up

1. Press the power button once, the alcohol tester LCD will turn on with one beep sound.

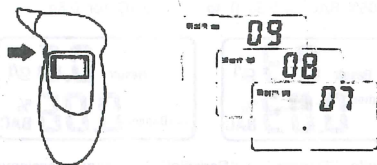2. The "warm up" symbol is displayed, and the 10 to 00 countdown is started that indicate the tester is in warm-up mode.

Note


* If first warm-up is defeated, the tester will be reset automatically and warm up again for 1-3 times.

* It always has a large number of gas & impurity adsorbed on the surface of alcohol sensor, when first time taking out the tester from the packaging, or long time no using the tester,

it need more time to clean and warm up the sensor, you should repress the power button and retry several times if possible.

Test

1. The "Blow" symbol is displayed when warm-up is succeed, and the 10 to 00 countdown is started that indicate the tester is in testing mode.

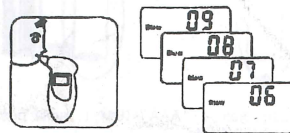2. Near and blow into the breath inhaler for 3-5 seconds.
3. Read the testing result from the alcohol tester LCD.
4. The "Caution" symbol is displayed if the alcohol concentration is the level:

- ☒ 0.02% BAC(or 0.2g/l) to 0.05% BAC(or 0.5g/l).
- ☐ 0.05% BAC(or 0.5g/l) to 0.08% BAC(or 0.8g/l).

5. The "Danger" symbol is displayed with beep audible alert if the alcohol concentration is the level:

- ☒ over & equal 0.05% BAC(or 0.5g/l).
- ☐ over & equal 0.08% BAC(or 0.8g/l).

Auto power off

1. After the result be displayed for about 10 seconds, the "OFF" symbol will appear for 2-3 seconds and the alcohol tester will be turned off.

2. If warm-up is defeated, the tester will be automatically turned off.
3. If no blow into the breath inhaler within 10 seconds in testing mode, the alcohol tester will be turned off after the LCD display the result "0.00% BAC & 0.0g/l" for about 10 seconds.

Note

* To obtain more accurate result, manufacturer recommend you to check repeatedly 2 or 3 times.

Replace battery

The LCD flickering indicates low power of batteries, please replace two "AAA" alkaline batteries.


Specifications

- * Sensor: advanced semiconductor oxide alcohol sensor
- * Range: 0.00%-0.19% BAC (0.0-1.9g/l)
- * Power input: 2x1.5V "AAA" alkaline battery
- * Operating time: warm up: 10s, response: <5s
- * Dimension: 95x36x18mm

Precaution

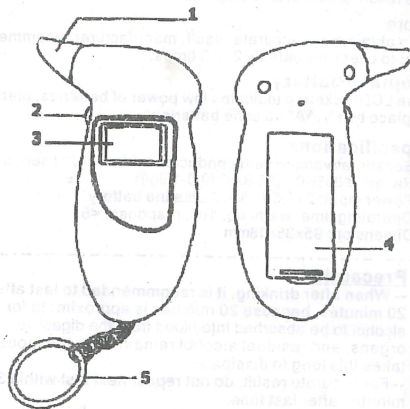- When after drinking, it is recommended to last after 20 minutes, because 20 minutes is approximate for alcohol to be absorbed into blood from the digestive organs, and residual alcohol remaining in the mouth takes this long to dissipate.
- For accurate result, do not repeat next test within 3 minutes after last time.
- Do not smoke and eat for 20 minutes prior testing.
- Do not use chemicals to clean the tester directly.
- Do not blowing liquid into the breath pipe.
- Do not use any excitant substances near the tester, such as paint, insecticides, alcohols.
- Do not place the tester in the closed environment with contaminated air.
- The result of this alcohol gives BAC for reference only, manufacturer does not to take any legal responsibility.

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